





RUCK LEG SUPPORT

RUCK®
FÜR FUSS UND PFLEGE

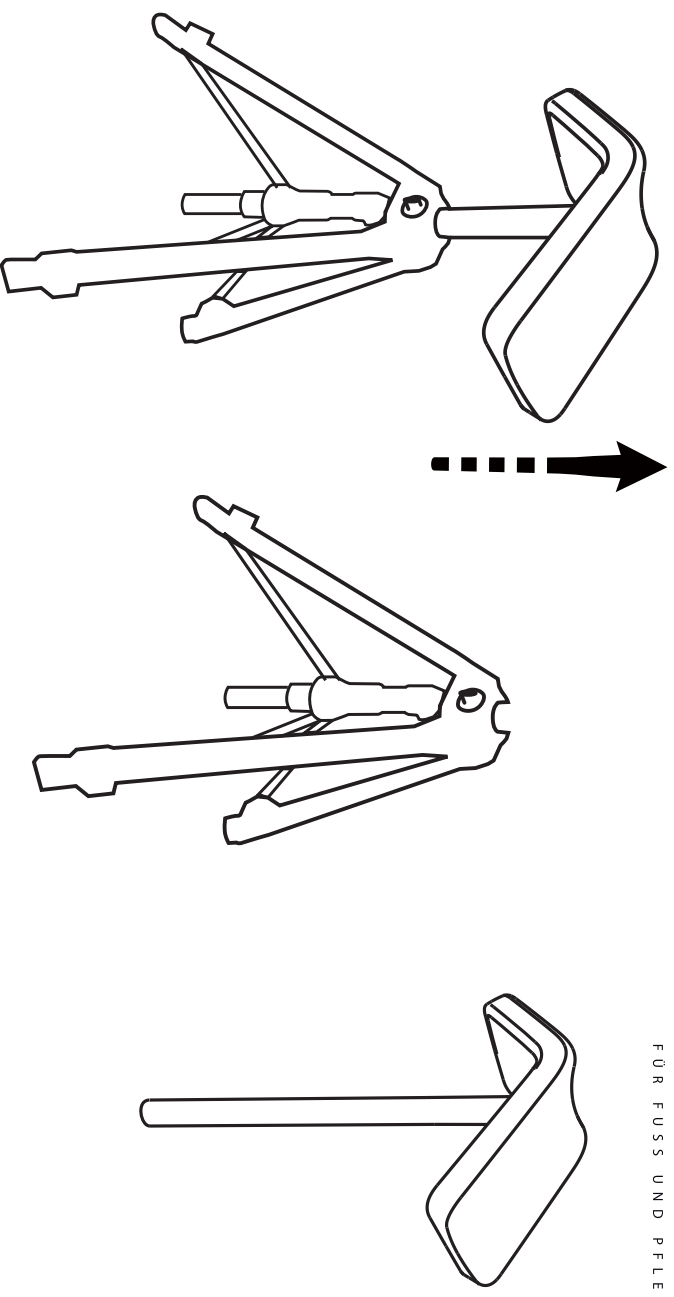


Operating instructions


Revision 02/2016 | Print

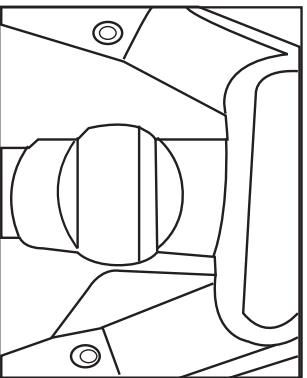


The quick adjustment lever allows for quick and easy height and position adaptation.

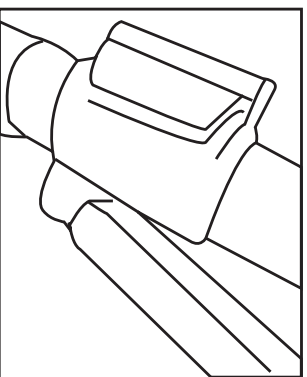


Use the lever to lock the tripod frame into position.

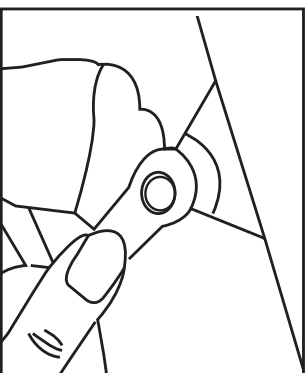

To open, move the clamping lever to the right, to tighten to the left.

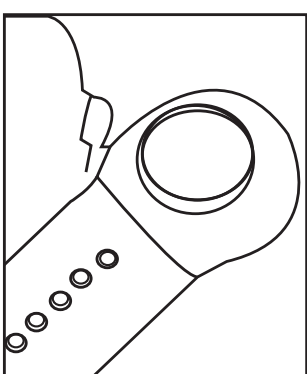

Only press the release button at the lever, if required.

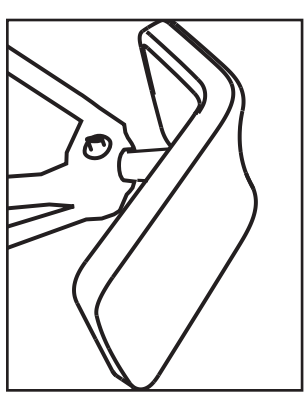

Support plate

Adjust the support plate using the clamping lever. To open the clamping lever move it to the right. For this, you do not need to press the release button at the centre of the lever. Use this button to adjust the position of the clamping lever. (For example, if the lever cannot be turned because it is too close to the plate).